## INSTRUCTIES

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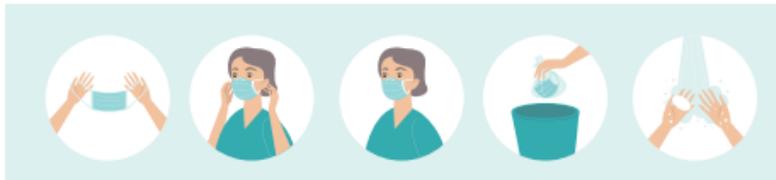
## INSTRUCTIONS

Description: community masks, pleated, non-sterile and meeting the requirements of standard NBN / DTD S 65-001: 2020 "Community masks and handmade masks - Guide to minimum requirements, making, maintenance and use"

Duration of use and precautions for use: single use for 4 hours maximum and throw away after use. If a mask is damaged or at the end of its life, it can, after having been washed, be thrown away with the usual waste. Unwashed masks should be wrapped in a closed plastic bag and then discarded in a closed bin. If a bin for medical waste is available, an unwashed used mask can be discarded.

Warnings: This product is not a medical device within the meaning of Regulation (EU) 2017/745 & Directive 93/42 / EEC or personal protective equipment within the meaning of Regulation (EU) 2016/425.

This product can only be used by children over 3 years of age due to the risk of suffocation and children should be supervised when wearing this product. Wearing a mask requires adaptation, namely the feeling of having less air, this is normal, the mask works. Please note: a mouth mask certainly does not allow intense effort, like intensive sport. This product is not designed to protect the lungs from hazards in the air. It is designed to avoid touching your mouth and nose and to minimize saliva splashes when talking, coughing or sneezing. The mask absolutely does not exempt the wearer from the application of barrier gestures supplemented by the measure of social distancing which are essential.

  
**MASKS FOR BELGIUM**  
[www.masksforbelgium.com](http://www.masksforbelgium.com)